

New!
One-of-a-Kind
3 in1 System



Great White® HD 1000 Stirrup Boot System

Introducing the World's First Interchangeable Stirrup Boot System with 1000 lb Patient Wt. Capacity

Instead of buying three sets of stirrups at a cost of \$15,000, buy our Great White® HD 1000 Interchangeable Stirrup Boot System and **Save 50%**!

Features

- **New** engineered design includes a stronger stainless steel tube for increased 1000 lb load capacity
- **New** redesigned 3-axis control handle mimics the anatomy of the human hip for smooth and accurate high/low lithotomy, adduction & abduction positioning
- **New** locking gas spring locks the stirrup in place and can not be moved unless the trigger is pulled
- Lithotomy Range -40° low to 75° high; Adduction -9° and abduction +25°
- Included Premier boot liners wrap around the foot, ankle & calf provide safe, comfortable & secure patient positioning

#800-0364 Great White® HD 1000 Stirrup Boot System

Required Accessory: Purchase Two SpringLoc Clamps

#800-0338 US #800-0338-EU EU
 #800-0338-JPN JPN #800-0338-SWISS SWISS


Optional Accessories

- #800-0074-S Great White® Stirrup Dolly
- #800-0284 Stirrup Wall Rack
- #800-0160 Clear View Leg Drapes
- #508-1415 Repl. Great White Boot Pads w/Fins
- #508-1354 Repl. Great White Boot Pads w/o Fins
- #508-0559 Repl. Great White Maxima Pads




System includes three pairs of Stirrup Boots:
 Platinum
 Premium
 Bariatric

Length & lithotomy indicators

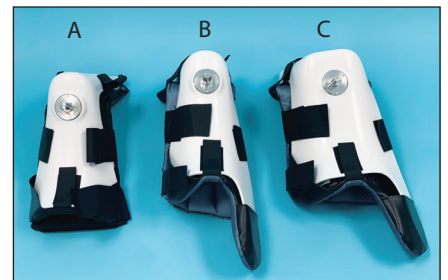

Select from the included Premium (A), Platinum (B) and Bariatric (C) Boots.

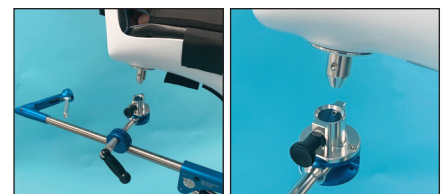

Quick and easy connection reduces setup time making it a cost effective product.

Phone 888.724.8763 • 781.982.7000 | Fax 781.982.7001 | www.schuremed.com

orders@schuremed.com | 452 Randolph Street, Abington, MA 02351