




Schure Grip Leg Holder

Provides versatile & excellent leg positioning

This surgical leg holder offers optimal patient positioning with dual independent braces that freely rotate 360° and spread from 3" to 11"W to fit any size leg, at any angle.

Features

- To use, twist one brace handle down, position the patient's leg on the leg holder, twist the handle back up, slide the brace horizontally for a snug fit, and lock in place with convenient twist handle
- Use also for ACL procedures
- Braces measure 12"L x 1.5"D x 8"H
- Mounting post measures 16"H x 5/8"D
- Required Schure Socket XL sold separately
- Stainless steel

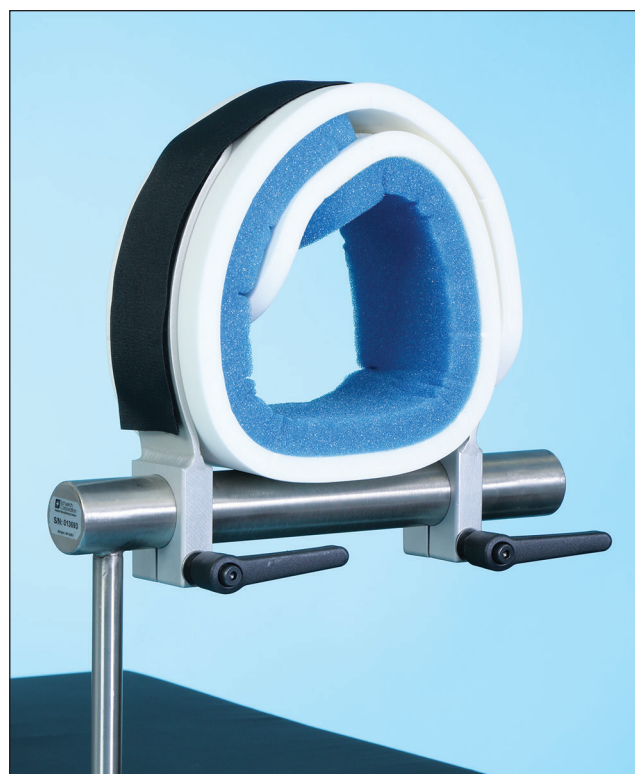

#800-0065 Schure Grip Leg Holder

Schure Socket XL

#800-0134 US

#800-0134-EU EU

#800-0134-JPN JPN

#800-0134-UK UK

#800-0134-DEN DEN

Optional Accessory

#508-0103 Disposable Premium Pads, 12/cs

Phone 888.724.8763 • 781.982.7000 | Fax 781.982.7001 | www.schuremed.com

orders@schuremed.com | 452 Randolph Street, Abington, MA 02351