

ROBOTIC STIRRUPS

#800-0342-RF (W/ FINS) | #800-0342-R (W/O FINS)



600 lb (272 kg)
Patient Weight Capacity

Lightly squeeze the ergonomic handle to easily allow full range of interoperative patient positioning



3-Axis control



Lower lithotomy range reduces the risk of positioning damage

Required Accessory
#800-0338
Available in
EU-JPN-SWISS



ROBOTIC STIRRUPS



Designed hand-in-hand with a robotic surgeon, SchureMed's Robotic Stirrups make robotic gynecologic & urologic procedures safer than they've ever been.

The mounting blade design quickly and easily inserts into side-rail clamps. The Great White Robotic Stirrups' range of motion is 60° low lithotomy to 60° high lithotomy. Adduction 9° and abduction 25° and also externally rotates out 20°, with the control to lock the medical stirrup for positioning anywhere in between. Two required SpringLoc Clamps sold separately.


Features

- ◆ Lithotomy indicators help to accurately reposition legs
- ◆ Simply squeeze the trigger for smooth and accurate high/low lithotomy, adduction and abduction positioning
- ◆ Lithotomy Range -60° low to 60° high; Adduction -9° and abduction +25°
- ◆ Included premier boot liners wrap the foot, ankle & calf providing safe, comfortable & secure patient positioning
- ◆ Exclusive 3-axis control from the handle mimics the anatomy of the human hip

ACCESSORIES



Wall
Rack
#800-0284

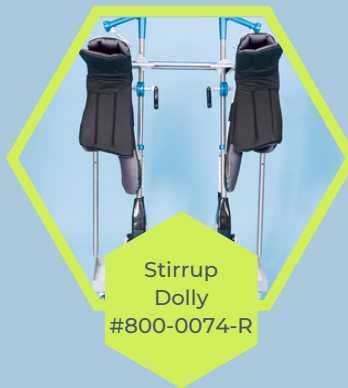

Stirrup
Dolly
#800-0074-R



Replacement
Pads
#508-1415

#800-0342-R/RF Great White Robotic Stirrups

Required Accessory:

Purchase Two SpringLoc Clamps

#800-0338 #800-0338-JPN
#800-0338-EU #800-0338-SWISS

Replacement Pads

#508-1415 Replacement Boot Pads w/ Fins
#508-1354 Replacement Boot Pads w/o Fins

Optional Accessories

#800-0074-R Robotic Stirrup Dolly
#800-0284 Great White Stirrup Wall Rack
#800-0160 Disp. Clear View Leg Drapes

For more information, or to place an order, contact our friendly and knowledgeable customer service team at +781-982-7000 or email us at orders@schuremed.com or visit us online at www.schuremed.com