



We put you and your patients in a better position



SURGICAL TABLE CLAMPS

888-724-8763 781-982-7000 | www.schuremed.com
452 Randolph Street, Abington, MA 02351



Established in 1994, SchureMed was founded by company President and CEO, Peter Schuerch. Peter's talent for problem solving led him from the operating room to the board room. Peter is not only the President and CEO, but also an integral part of the design team.

SchureMed manufactures high quality surgical assistance devices for medical facilities all over the world. Up or down, big or small, SchureMed is putting you, and your patients, in a better position.



OUR MISSION

Our mission is to enhance the operating room experience of patients and medical professionals worldwide. Our passion for innovation, quality, and safety drive us to create highly functional surgical positioning solutions that optimize results and reduce operating time through greater efficiency.

OUR VISION

Our vision is to improve quality of life for all patients and medical professionals by becoming the world's most trusted leader in innovative surgical positioning.



TABLE OF CONTENTS

02	Clamp List	03
Schure Socket XL		
04	Schure Socket	05
Stirrup Clamp		
06	SpringLoc Clamp	07
CamLoc Rail Clamp		
08	Deluxe Rail Clamp	09
Universal Rail Clamp		
10	Simple Rail Clamp	11
Infinity Line		

Schure Socket XL

800-0134
1000 lbs (453 kg)
1" x 3/8" (2.5cm x .95 cm)
3/4"D (1.9 cm)




Schure Socket

800-0006
800 lbs (363 kg)
5/8"D (1.6 cm)



CamLoc Clamp

800-0271
500 lbs (227 kg)
5/8"D (1.6 cm)
1" x 3/8" (2.5 cm x .95 cm)

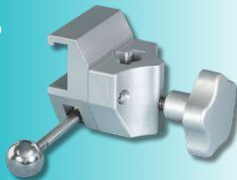

Deluxe Clamp

800-0248
500 lbs (227 kg)
5/8"D (1.6 cm)
1" x 3/8" (2.5 cm x .95 cm)



Universal Clamp

800-0085
350 lbs (159 kg)
5/8"D (1.6 cm)
1" x 3/8" (2.5 cm x .95 cm)



Simple Clamp

800-0228
720 lbs (327 kg)
1" x 3/8" (2.5cm x .95cm) vertical
1" x 1/4" (2.5cm x .635cm) horizontal



Stirrup Clamp

800-0396
1000 lbs (453 kg)
1" x 3/8" (2.5 cm x .95 cm)



SpringLoc Clamp

800-0127
800 lbs (363 kg)
1" x 1/4" (2.5 cm x .635 cm)



Infinity

800-0386
800 lbs (363 kg)
5/8"D to 3/4"D (1.6 cm to 1.9 cm)
1" x up to 3/8" (2.5 cm x .95 cm)



Infinity sq.

800-0386-SQ
800 lbs (363 kg)
.787" (20 mm) square bar



Infinity Fixed

800-0386-F
800 lbs (363 kg)
5/8"D to 3/4"D (1.6 cm to 1.9 cm)
1" x up to 3/8" (2.5 cm x .95 cm)

The revolutionary Schure Socket XL snaps on at any rail location and fits all standard operating room table side rails measuring 3/8"W x 1-1/8"H (.95 cm x 3 cm)

Hardened clamping surface creates maximum clamping force, stopping all movement of operating room accessories

Accepts rectangular mounting posts up to 1" x 3/8" (2.5 cm x .95 cm)

Accepts any round mounting posts up to 5/8" (1.6 cm) in diameter

Patient weight capacity is 800 lbs (363 kg)

Quickly and easily connects and disconnects securely anywhere along the side rail with a quick turn, even over side rail stand off

NEXT
Generation



Schure Socket XL

#800-0134
#800-0134-DEN

#800-0134-EU
#800-0134-JPN

#800-0134-UK
#800-0134-SWISS

For use with:



Operating room staff can quickly position surgical devices anywhere along standard operating room table side rails - Use the Schure Socket to prevent all movement on mounting posts



Fits round mounting posts up to 5/8" (1.6 cm) in diameter

Patient weight capacity is 500 lbs (227 kg)

The industrial-strength Schure Socket fits all operating room tables with standard 3/8"W x 1-1/8"H (.95 cm x 2.85 cm) side rails

Connects anywhere along side rail, without the need for notches



Schure Socket

#800-0006
#800-0006-DEN

#800-0006-EU
#800-0006-JPN

#800-0006-UK
#800-0006-SWISS

For use with:



Our low-profile Stirrup Clamp is one of the world's easiest-to-use blade clamps, effortlessly securing patient-positioning devices and surgical equipment with a flat mounting post anywhere along an operating-table side rail

Securely locks accessories with a flat 1" x 1/4" (2.5 cm x .64 cm) mounting post

1,000 lbs (454 kg) patient weight capacity

Drops effortlessly onto the operating table side rail without the need for table notches

Easily removeable without the need for a release lever




Stirrup Clamp

#800-0396
#800-0396-DEN

#800-0396-EU
#800-0396-JPN

#800-0396-UK
#800-0396-SWISS

For use with:

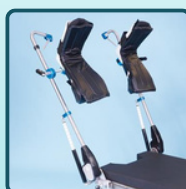

800 lbs (363 kg)
patient weight
capacity

Springloc Clamp is easily
removed with gentle
pressure applied to the lever



Securely locks
accessories
with a flat 1" x 1/4"
(2.5 cm x .64 cm)
mounting post

Drops effortlessly onto
the operating table side
rail without the need for
table notches


NEXT
Generation



SpringLoc Clamp

#800-0127


For use with:

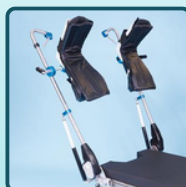

Accepts accessories with round posts up to 5/8" (1.6 cm) in diameter and flat mounting posts to 1" x 3/8" (2.5 cm x .95 cm)

Attaches to the operating table rail at any location

Holds 1-2 drapes between table side rail and clamp

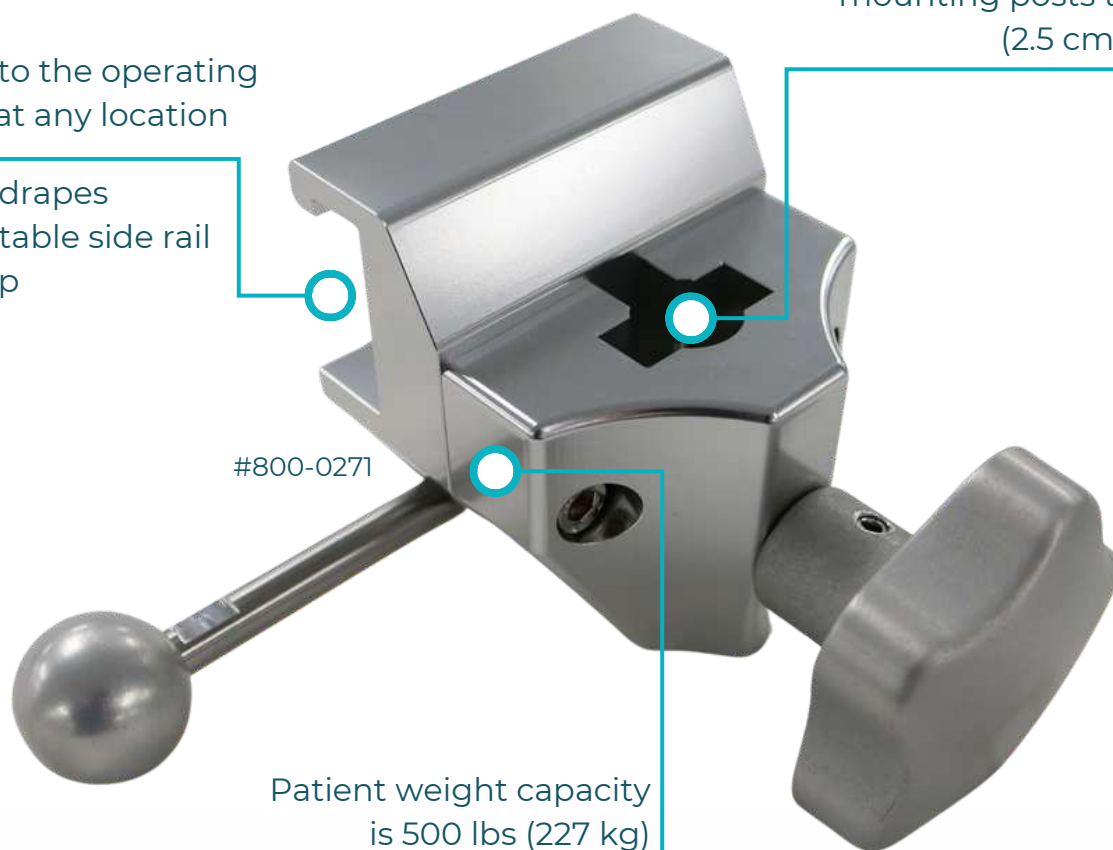

Available in low profile - a more sleek and compact design
#800-0271-LP



Available in large - for more drape capacity (Holds 3-5 drapes between table side rail and clamp)
#800-0271-L

Designed for use in the sterile field for when drapes are to be placed between the clamp and siderail to maintain sterility; attaching common accessories to nearly all OR table designs and brands



CamLoc Rail Clamp

#800-0271-DEN
#800-0271-EU

#800-0271-JPN
#800-0271-UK

#800-0271-SWISS
#800-0271-LP-DEN

#800-0271-LP-EU
#800-0271-LP-JPN

#800-0271-LP-UK
#800-0271-LP-SWISS

For use with:





Designed for a wide-range of purposes; affixing most operating room table accessories to most OR table designs and brands

Accessories can be mounted either vertically or horizontally

Securely locks all accessories with 1" x 3/8" (2.5 cm x .95 cm) flat mounting post

Patient weight capacity is 500 lbs (227 kg)

Attaches at OR table rail ends or notches

Securely locks all accessories with round mounting posts up to 5/8" (1.6 cm) in diameter




Deluxe Rail Clamp

#800-0248
#800-0248-DEN

#800-0248-EU
#800-0248-JPN

#800-0248-UK
#800-0248-SWISS

For use with:

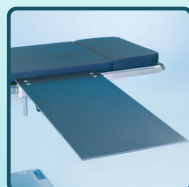

Designed for a wide-range of purposes; affixing most operating room table accessories to most OR table designs and brands

Accepts any round mounting posts up to 5/8" (1.6 cm) in diameter and rectangular mounting posts up to 1" x 3/8" (2.5 cm x .95 cm)



Patient weight capacity is 350 lbs (159 kg)

Snaps on at any rail location and attaches to all 3/8"W x 1-1/8"H (.95 cm x 3 cm) surgical table side rails

Holds accessories in the vertical and horizontal position



Universal Rail Clamp


#800-0085

For use with:



Designed to secure accessories with a flat mounting blade to most OR table designs and brands

Patient weight capacity is 720 lbs (327 kg)

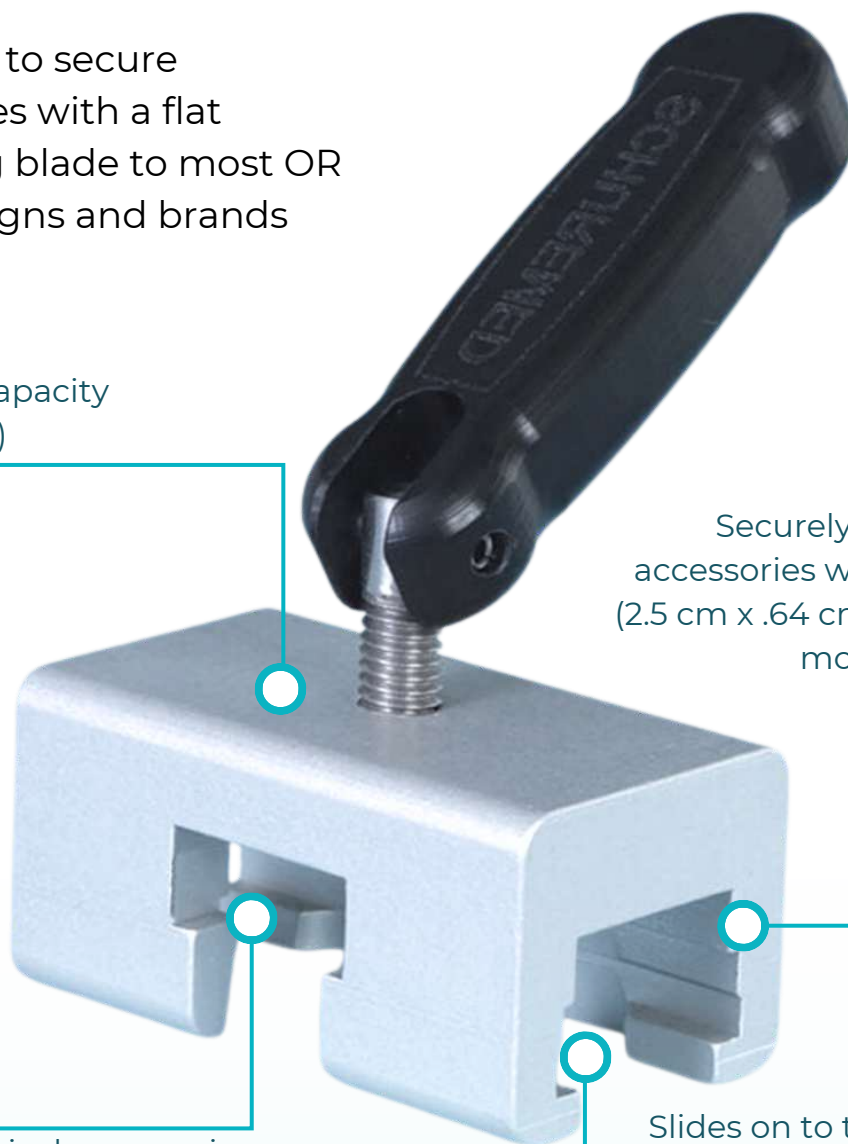

Securely locks surgical accessories with flat 1" x 1/4" (2.5 cm x .64 cm) horizontally mounting blades

Securely locks surgical accessories with flat 1" x 3/8" (2.5 cm x .95 cm) vertically mounting blades

Slides on to the rail at operating table rail end or notches






Simple Clamp

#800-0228
#800-0228-DEN

#800-0228-EU
#800-0228-JPN

#800-0228-UK
#800-0228-SWISS

For use with:

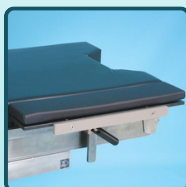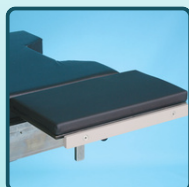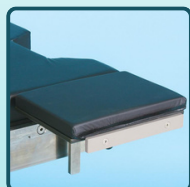



The Infinity Line's adjustable jaw allows Infinity Clamps to attach to any OR table in the world – eliminating the need to purchase different clamps for different surgical table side rail sizes



Infinity

#800-0386

5/8" (16 mm) Diameter Round Rod
.709" (18 mm) Diameter Round Rod
3/4" (19 mm) Diameter Round Rod
3/8" x 1" (.95 cm x 2.5 cm) Square Rod




Infinity sq.

#800-0386-SQ

.787" (20 mm) Square Rod

Infinity Line Clamps have
a 800 lbs (363 kg)
Patient Weight Capacity

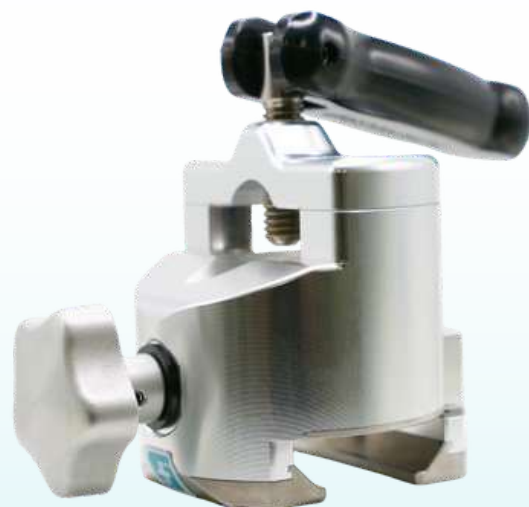

Infinity Fixed

#800-0386-F


Stirrup Blade
5/8"D (16 mm) Round Rod
.709" (18 mm) Round Rod
3/4"D (19 mm) Round Rod
1/4" x 1" (.635 cm x 2.5 cm) Flat Bar
3/8" x 1" (.95 cm x 2.5 cm) Flat Bar



Infinity Line Clamps

#800-0386
#800-0386-F

#800-0386-SQ





+1-888-724-8763



orders@schuremed.com



452 Randolph Street
Abington, MA, USA 02351



www.schuremed.com



Rev. 001.24