

GENOS TKR POSITIONER

#800-0409 30" GENOS TKR Positioner

Pull apart system
allows for smaller
sterilization tray

Place the operative leg in any
degree of flexion to full
extension

The ball track system at the
heel of the leg holder allows
surgeons to tilt and rotate the
operative leg

Track locking bracket securely
locks the Base Plate to the
table with two thumb screws

Built-in stops on each end of
the track system prevent the
leg holder from sliding off of
the Base Plate

Locks in any position
with a simple turn of
the wrist

Smooth sliding motion for
effortless leg positioning




INCLUDED ACCESSORIES



All components are autoclavable (Exposure Time at 270°F (132°C) for 4 minutes)



#800-0409-TRAY
Sterilization Tray

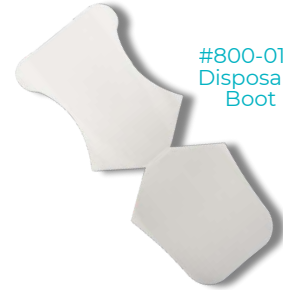

#800-0146-SGL
Disposable TKR
Boot Liner

Included:
1x TKR Positioner
1x TKR Disposable Boot Liner
1x TKR Sterilization Tray

REQUIRED ACCESSORIES



#800-0271
CamLoc Clamp

Required Accessories
#800-0271 CamLoc Clamp

OPTIONAL ACCESSORIES

Disposables
#800-0146 12x Disposable TKR Boot Liners



#800-0146
Case of 12
Disposable TKR
Boot Liners

Optional Accessories
#800-0271-L CamLoc Clamp Large
#800-0271-LP CamLoc Clamp Low-Profile



#800-0271-L
CamLoc Clamp
Large



#800-0271-LP
CamLoc Clamp
Low-Profile